CE Diameter Report - Number Distribution



Measurement Details

Sample Name: FILLITE EGP5

SOP Name: MATERIE PRIME.vsop **Edited:** False

User Name: morphologi Edited On: 06 May 2014 10:22:51

Date: 06 May 2014 10:22:51

Analysis

SOP Analysis: 3.0 SOP Fill Holes: True

SOP Segmentation: None **SOP Trash Size**: 10

SOP Optic(s) used: 2.5x Particles Counted: 286317


Result

CE Diameter Minimum (µm): 4.29 CE Diameter D[n, 0.1] (µm): 4.43

CE Diameter Maximum (µm): 649.25 CE Diameter D[n, 0.5] (µm): 6.85

CE Diameter Mean (µm): 8.39 CE Diameter D[n, 0.9] (µm): 13.16

CE Diameter STDV (µm): 10.76 CE Diameter RSD (%): 128.13

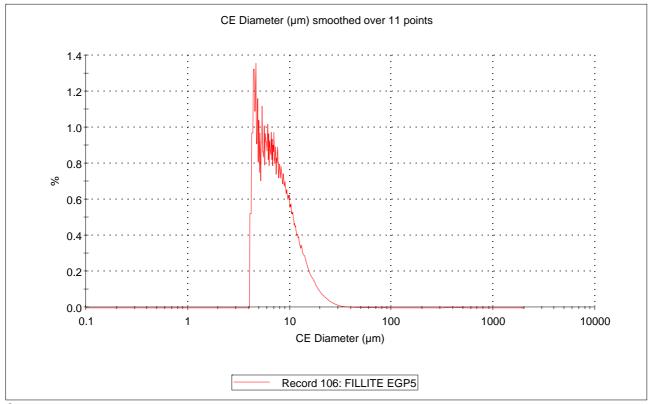

Comments:

Morphologi Malvern Instruments Worldwide www.malvern.com