CE Diameter Report - Number Distribution



Measurement Details

Sample Name: Italcementi grigio 42.5

SOP Name: 20X.vsop **Edited:** False

User Name: Morphologi Edited On: 21 October 2015 11:20:12

Date: 21 October 2015 11:20:12

Analysis

SOP Analysis: 3.0 SOP Fill Holes: True

SOP Segmentation: None **SOP Trash Size**: 10

SOP Optic(s) used: 20x Particles Counted: 451551


Result

CE Diameter Minimum (µm): 0.54 CE Diameter D[n, 0.1] (µm): 1.31

CE Diameter Maximum (µm): 120.46 CE Diameter D[n, 0.5] (µm): 3.33

CE Diameter Mean (µm): 4.33 CE Diameter D[n, 0.9] (µm): 7.95

CE Diameter STDV (µm): 3.87 CE Diameter RSD (%): 89.45

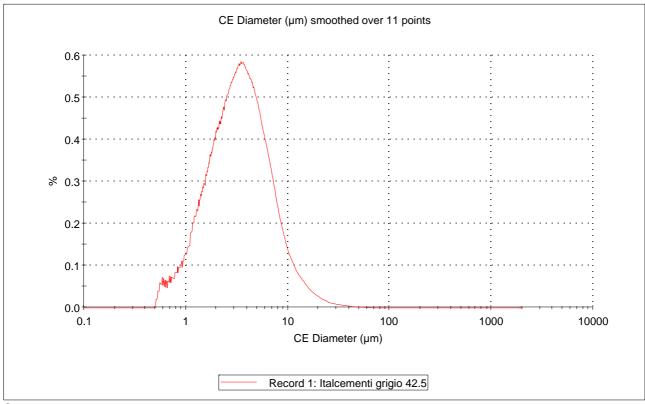

Comments:

Morphologi Malvern Instruments Worldwide www.malvern.com