

CE Diameter Report - Volume Distribution



Measurement Details

Sample Name: FILLITE EGP5

SOP Name: MATERIE PRIME.vsop

Edited: False

User Name: morphologi

Edited On: 06 May 2014 10:22:51

Date: 06 May 2014 10:22:51

Analysis

SOP Analysis: 3.0

SOP Fill Holes: True

SOP Segmentation: None

SOP Trash Size: 10

SOP Optic(s) used: 2.5x

Result

CE Diameter Minimum (μm): 4.29

CE Diameter D[v, 0.1]: 358.6

CE Diameter Maximum (μm): 649.25

CE Diameter D[v, 0.5]: 429.9

CE Diameter D[4,3] (μm): 438.5

CE Diameter D[v, 0.9]: 538.1


CE Diameter D[3,2] (μm): 244.5

CE Diameter STDV (μm): 10.76

Particles Counted: 286317

CE Diameter RSD (%): 128.13

Volume transformation: CE Diameter (μm) smoothed over 11 points

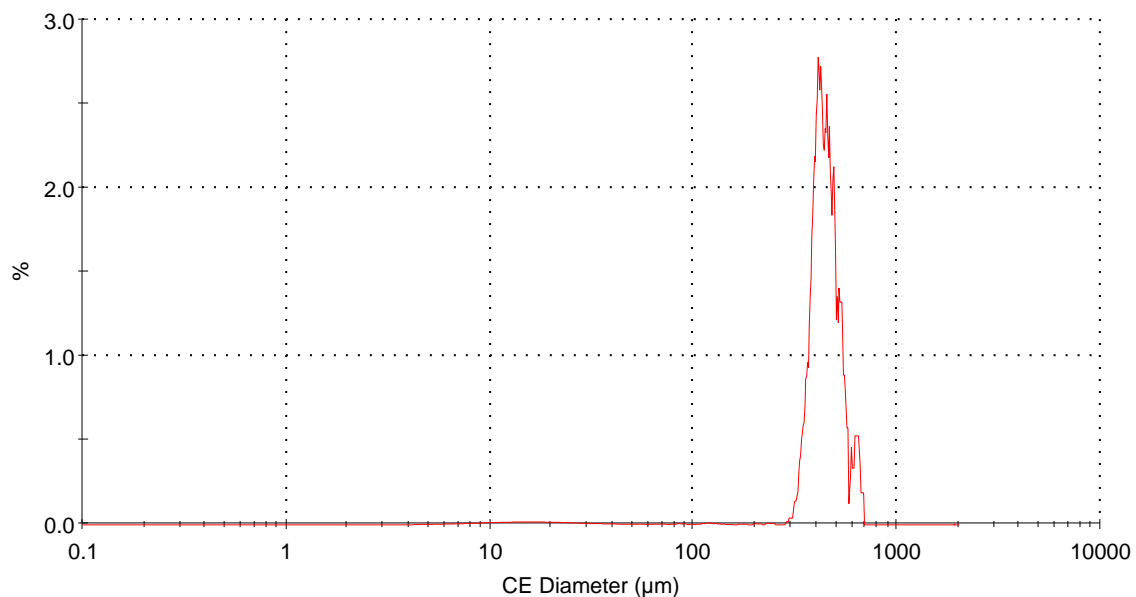

Record 106: FILLITE EGP5

Comments:

Morphologi

Malvern Instruments Worldwide

www.malvern.com