

CE Diameter Report - Number Distribution



Measurement Details

Sample Name: DRY D1 C NEW GENERATION

SOP Name: MATERIE PRIME.vsop

Edited: True

User Name: Morphologi

Edited On: 29 April 2014 08:30:19

Date: 28 April 2014 20:09:44

Analysis

SOP Analysis: 3.0

SOP Fill Holes: True

SOP Segmentation: None

SOP Trash Size: 10

SOP Optic(s) used: 20x, 50x

Particles Counted: 19058

Result

CE Diameter Minimum (μm): 0.22

CE Diameter D[n, 0.1] (μm): 1.49

CE Diameter Maximum (μm): 54.31

CE Diameter D[n, 0.5] (μm): 3.76


CE Diameter Mean (μm): 4.45

CE Diameter D[n, 0.9] (μm): 7.70

CE Diameter STDV (μm): 3.31

CE Diameter RSD (%): 74.32

CE Diameter (μm) smoothed over 11 points

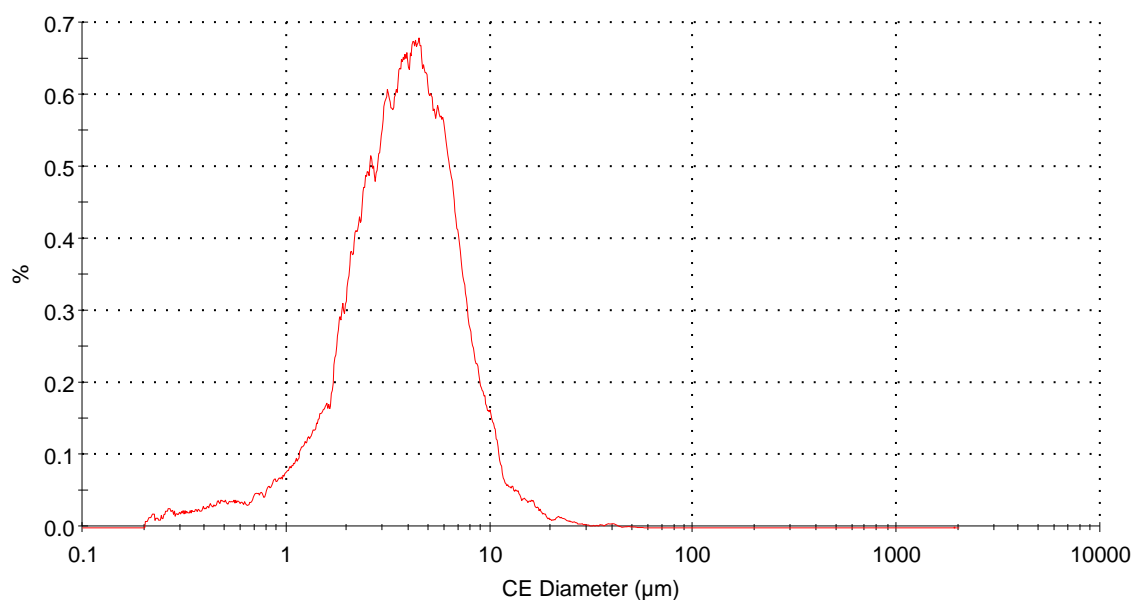

Record 82: DRY D1 C NEW GENERATION

Comments:

Morphologi

Malvern Instruments Worldwide

www.malvern.com