

CE Diameter Report - Volume Distribution



Measurement Details

Sample Name: Filler Marmifera mineraria

SOP Name: MATERIE PRIME.vsop

Edited: False

User Name: Morphologi

Edited On: 27 April 2015 11:55:51

Date: 27 April 2015 11:55:51

Analysis

SOP Analysis: 3.0

SOP Fill Holes: True

SOP Segmentation: None

SOP Trash Size: 10

SOP Optic(s) used: 50x

Result

CE Diameter Minimum (μm): 0.22

CE Diameter D[v, 0.1]: 6.834

CE Diameter Maximum (μm): 77.97

CE Diameter D[v, 0.5]: 15.62

CE Diameter D[4,3] (μm): 21.06

CE Diameter D[v, 0.9]: 43.31


CE Diameter D[3,2] (μm): 12.95

CE Diameter STDV (μm): 4.13

Particles Counted: 27309

CE Diameter RSD (%): 89.26

Volume transformation: CE Diameter (μm) smoothed over 11 points

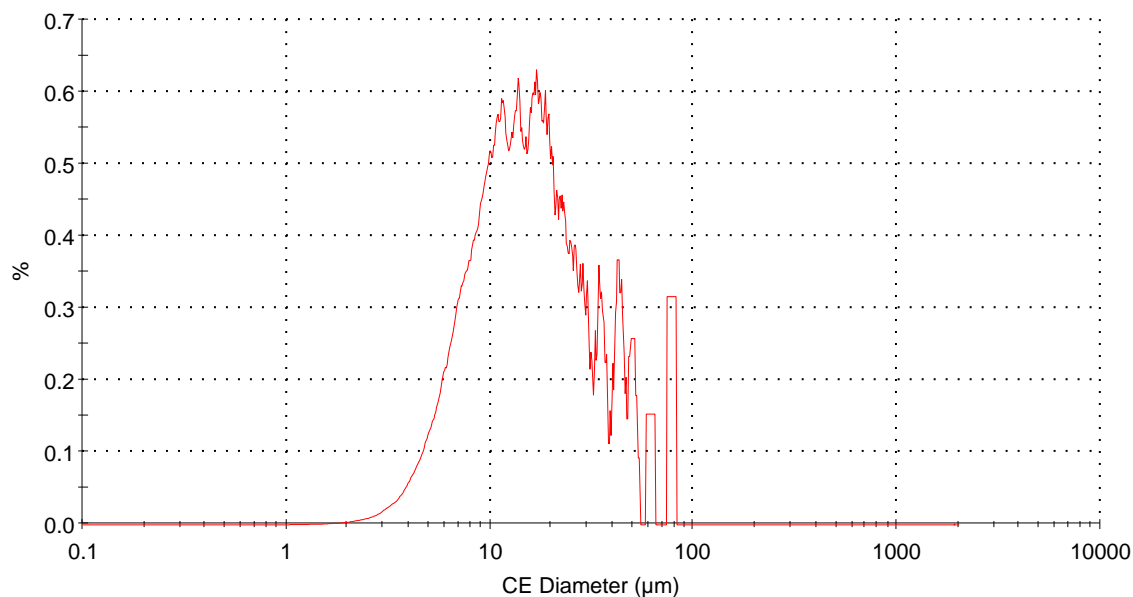

Record 19: Filler Marmifera mineraria

Comments:

Morphologi

Malvern Instruments Worldwide

www.malvern.com