

# CE Diameter Report - Volume Distribution



## Measurement Details

**Sample Name:** Formiato Shandong

**SOP Name:** 10X.vsop

**Edited:** False

**User Name:** Morphologi

**Edited On:** 27 November 2015 13:10:51

**Date:** 27 November 2015 13:10:51

## Analysis

**SOP Analysis:** 3.0

**SOP Fill Holes:** True

**SOP Segmentation:** None

**SOP Trash Size:** 10

**SOP Optic(s) used:** 10x

## Result

**CE Diameter Minimum ( $\mu\text{m}$ ):** 1.08

**CE Diameter D[v, 0.1]:** 54.27

**CE Diameter Maximum ( $\mu\text{m}$ ):** 321.95

**CE Diameter D[v, 0.5]:** 103.8

**CE Diameter D[4,3] ( $\mu\text{m}$ ):** 110.6

**CE Diameter D[v, 0.9]:** 169.5


**CE Diameter D[3,2] ( $\mu\text{m}$ ):** 88.49

**CE Diameter STDV ( $\mu\text{m}$ ):** 28.73

**Particles Counted:** 22677

**CE Diameter RSD (%):** 129.17

Volume transformation: CE Diameter ( $\mu\text{m}$ ) smoothed over 11 points

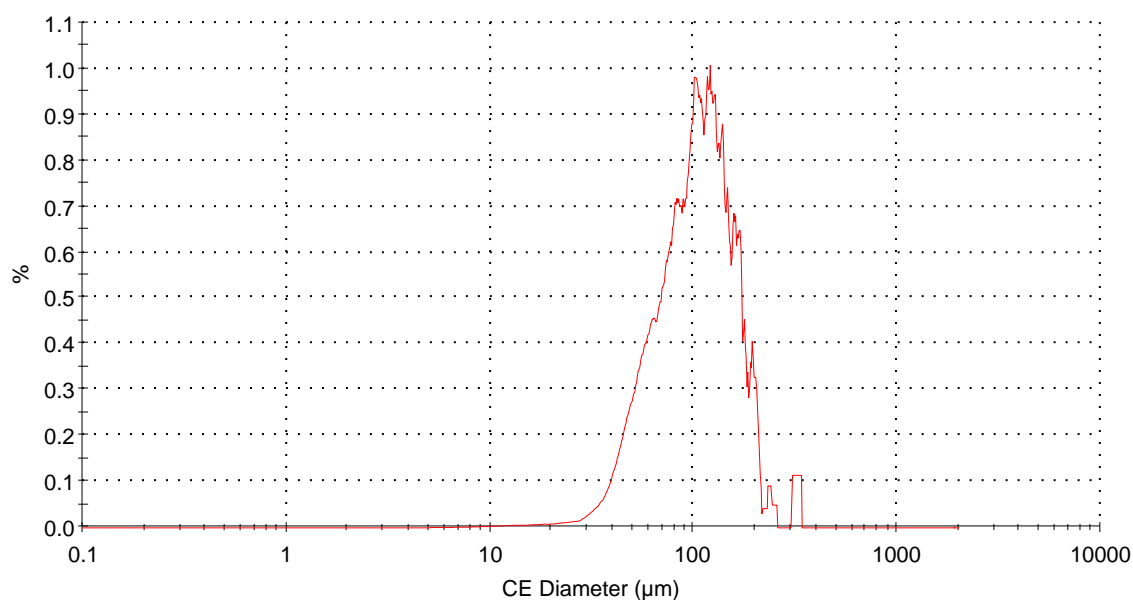

Record 1: Formiato Shandong

## Comments:

Morphologi

Malvern Instruments Worldwide

[www.malvern.com](http://www.malvern.com)