CE Diameter Report - Number Distribution



Measurement Details

Sample Name: ARBOCEL PWC 50-33 LD

SOP Name: MATERIE PRIME.vsop **Edited:** True

User Name: Edited On: 21 May 2014 18:49:37

Date: 21 May 2014 18:45:49

Analysis

SOP Analysis: 3.0 SOP Fill Holes: True

SOP Segmentation: None SOP Trash Size: 10

SOP Optic(s) used: 2.5x Particles Counted: 31938


Result

CE Diameter Minimum (µm): 4.29 CE Diameter D[n, 0.1] (µm): 4.61

CE Diameter Maximum (µm): 453.77 CE Diameter D[n, 0.5] (µm): 8.46

CE Diameter Mean (μm): 14.17 CE Diameter D[n, 0.9] (μm): 19.68

CE Diameter STDV (µm): 23.83 CE Diameter RSD (%): 168.13

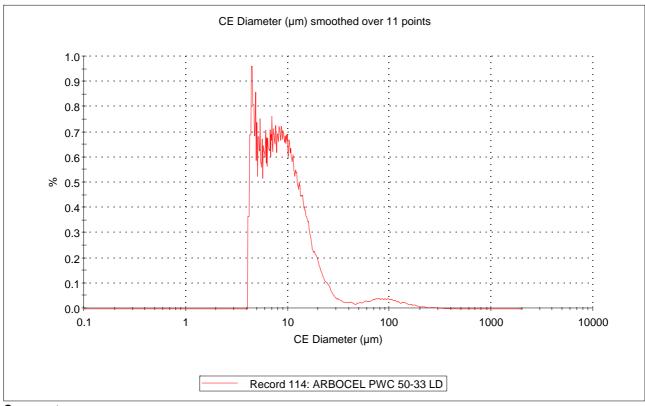

Comments:

Morphologi Malvern Instruments Worldwide www.malvern.com