CE Diameter Report - Number Distribution



Measurement Details

Sample Name: VAEPOL DMC CIMENTAL

SOP Name: MATERIE PRIME.vsop Edited: True

User Name: Edited On: 13 May 2014 08:41:08

Date: 13 May 2014 05:48:24

Analysis

SOP Analysis: 3.0 SOP Fill Holes: True

SOP Segmentation: None **SOP Trash Size**: 10

SOP Optic(s) used: 10x, 20x Particles Counted: 157906


Result

CE Diameter Minimum (µm): 0.54 CE Diameter D[n, 0.1] (µm): 1.17

CE Diameter Maximum (µm): 184.31 CE Diameter D[n, 0.5] (µm): 3.97

CE Diameter Mean (µm): 6.27 CE Diameter D[n, 0.9] (µm): 12.36

CE Diameter STDV (µm): 8.32 CE Diameter RSD (%): 132.83

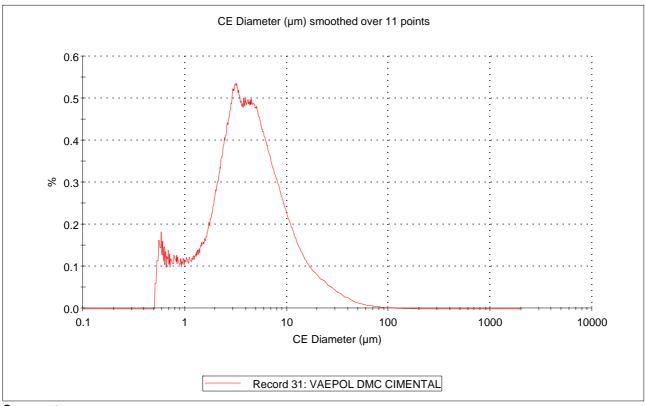

Comments:

Morphologi Malvern Instruments Worldwide www.malvern.com